



Locking the flow rate (Millennium concentrator)

Your child's concentrator will be locked at the prescribed flow rate by your Air Products technician. However, if your child has been prescribed a variable flow rate, you will need to adjust the flow rate and relock the machine. The following instructions explain how to do this.

How to lock the flow rate

- 1.** First locate the grub screw behind the flowmeter valve (A).
- 2.** Open the filter housing on the rear of concentrator.
- 3.** Remove the Allen key located in the sponge padding.
- 4.** Insert the Allen key into grub screw beneath the flowmeter valve (A).
- 5.** If the flow is already locked, unlock it by turning the Allen key anti-clockwise.
- 6.** Set the flow to desired flow rate.
- 7.** Re-lock the flow rate by turning the grub screw clockwise with the Allen key.
- 8.** Ensure the ball stays on the correct setting by adjusting the flowmeter valve (A) while tightening the grub screw.

Fig.B Locking the flow rate on the Millennium concentrator

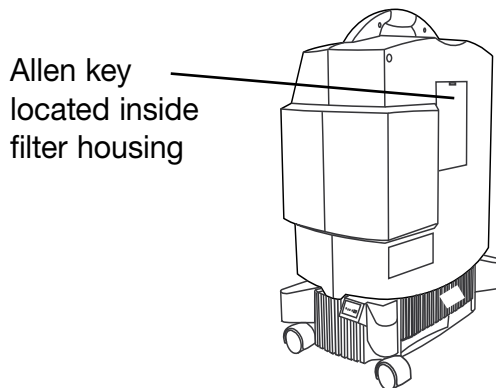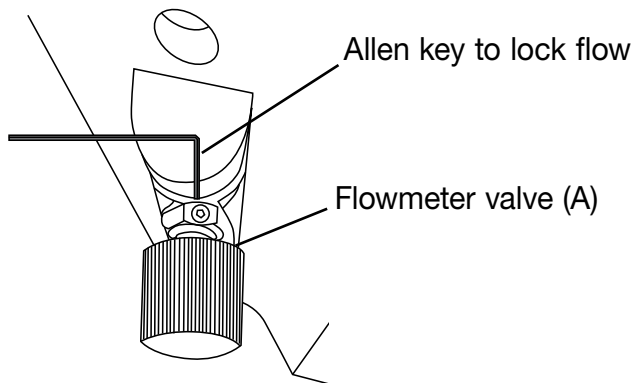

Equipment for Children and Babies